

Medium universal catheter securement with wide silicone adhesive area Grip-Lok®

Intended for use in the general population, for adhesive securement of the catheter hub, tube or line to the patient's skin.

1. **Select** the area for the placement of the Grip-Lok. Note: catheter, tube or line should be already placed/inserted.
2. **Prepare** the skin according to the standard hospital protocol for dressing application. Skin prep or hair removal may be required on some patients for better adhesion.
3. **Open** the top flap and slide the Grip-Lok, under the catheter, tube or line, centering it above the adhesive area.
4. **Pull** the paper backing from one side of the Grip-Lok, then the other, to secure in desired position on the skin.
5. **Remove** the interior liner to expose the adhesive areas.
6. **Place** the catheter, tube or line in the center of the exposed adhesive.
7. **Secure** top flap over the catheter, tube or line and apply gentle pressure to the center fabric area over the adhesive.

Direction: For additional skin protection, use skin prep pad prior to application of stabilization device.

Note: Use of an alcohol swab may help removal of securement device from the skin.

⚠ If you experience severe redness, itching, swelling, or irritation of the skin, consult your physician as this may be a sign of an allergic reaction


⚠ Re-use of this device may change its mechanical or biological features and may cause device failure, allergic reactions or bacterial infections.

⚠ Replace securement device if soiled or saturated in fluid or if device shows signs of wear or damage.

Note: if being used in a medical facility, dispose according to medical facility and municipal guidelines.

⊘ Do not use if the product's sterile barrier system or its packaging is compromised.

Report any serious incident related to device use to TIDI Products and the Member State competent authority.

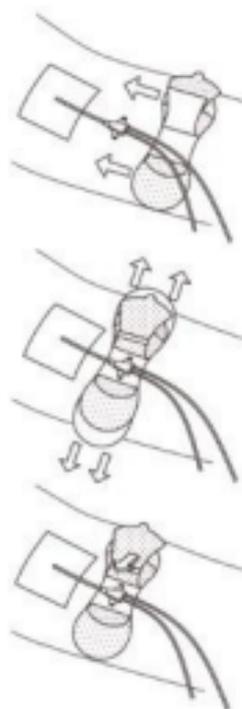

44U00010 Rev. 03 2023-08-21

HYPOALLERGENIC



TIDI® Products, LLC
570 Enterprise Drive
Neenah, WI 54956 USA

EC REP

MDSS GmbH
Schiffgraben 41
30175 Hannover · DE

STERILE R

CE 2797